

Quantifying T cells or natural killer cells secreting IFN- γ , TNF- α , or both cytokines in a multiplex assay using the Xdrop[®] workflow

Mette Juul Jacobsen, Mary Amasia, Christian Hussing, Rikke Stick Højmark, Thomas Mørch Frøsig, and Marie Just Mikkelsen
Samplix ApS, Birkerød, Denmark

Summary

- This Xdrop workflow enables rapid, multiplex quantification of the immune cells in a population that are secreting one or more cytokines.
- The multiplex and singleplex assay results concurred and were highly reproducible.

Introduction

Bulk cell analyses can miss crucial information on the homogeneity of immune cell populations and the diversity of individual immune cell responses. Within cell therapy research, there is a clear need for rapid, sensitive, and reliable functional analyses of cells with single-cell resolution to gain these critical insights.

Samplix has developed Xdrop and the Xdrop DE50 Cartridge to encapsulate living mammalian cells in highly stable double-emulsion droplets (DE50 droplets) for incubation, flow cytometry, and sorting.

The Xdrop workflow provides the required single-cell view of cytokine levels, because the droplet encapsulation prevents crosstalk between the cells. It allows for the analysis of 150,000 single cells per sample with the option to recover immune cells with the desired functionality for later molecular profiling or further cell expansion.

Furthermore, results are obtained in less than eight hours, which is compatible with work processes at most cell biology laboratories. This is thanks to the small droplet volumes (~100 pL), which ensure secreted cytokines rapidly reach detectable levels.

Here, we use the Xdrop workflow to encapsulate immune cells in droplets for the quantification of TNF- α and IFN- γ secretion and digital counting of the number of actively secreting cells. We also show that the workflow results are highly reproducible.

Experimental setup

Human natural killer cells (NK-92 cell line) and a population of human peripheral blood mononuclear cells (PBMCs) from a healthy donor were analyzed in parallel.

The NK-92 cells and PBMCs were labeled with IFN- γ and TNF- α capture reagents according to the Miltenyi Biotec[®] IFN- γ and TNF- α Secretion Assay protocols. The cell populations were divided in two: a non-activated control subpopulation that was not exposed to cytokine stimulants; and an activated population that was exposed to 25 ng/ml cytokine stimulation cocktail, 100 ng/ml IL-2, and 2.5 μ g/ml PHA-M. This was earlier shown to activate the NK-92 cells and the T cells among the PMBCs (data not shown).

The control and activated subpopulations of each cell type were then divided into samples of ~150,000 cells for encapsulation in DE50 droplets together with IFN- γ APC and TNF- α PE antibodies.

The cells within the DE droplets were subsequently incubated for 3 hours at 37°C in 5% CO₂. After incubation, all the cells were recovered from the droplets, washed, and stained with LIVE/DEAD Fixable Green Dead Stain (Thermo Fisher Scientific[®]) to identify living cells. The PBMCs were also stained with CD3-PerCp to identify the T cells within the population. All the cells were then analyzed using a BD Accuri[™] flow cytometer.

To verify that the multiplex assay yields comparable results to a singleplex assay, the experiment was repeated for each cell type and each cytokine singly, i.e., with one capture reagent bound to the cells instead of two.

The reproducibility of this Xdrop workflow for cytokine secretion analysis was determined by repeating the same experiments 5 times on 3 different days.

Figure 1 shows the basic workflow.

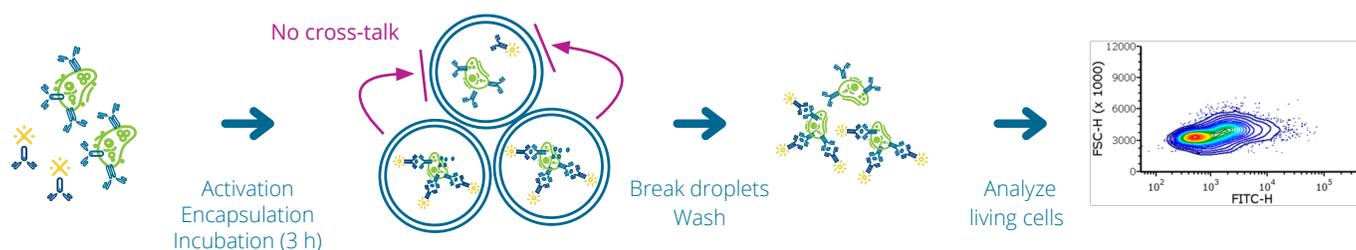


Figure 1. The Xdrop workflow for a single-cell cytokine secretion assay. The whole workflow takes ~6 hours based on the 3 hours of incubation used here. Two to four hours of incubation has been found to be sufficient in other experiments (data not shown).

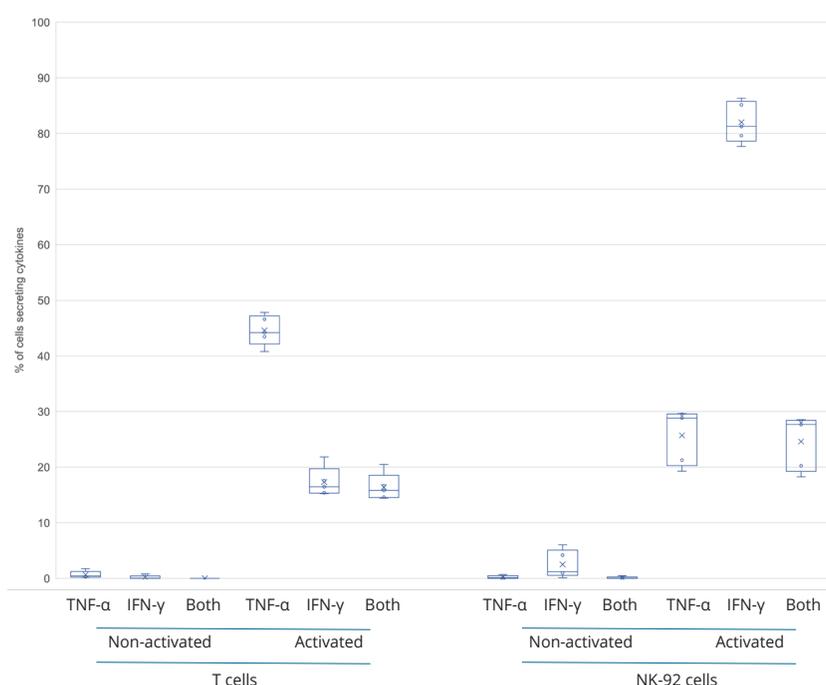


Figure 2. Box plot showing the TNF- α and IFN- γ secretion profiles for non-activated (control) and activated T cells within a PMBC population from a human blood sample and for non-activated (control) and activated NK-92 cells. Data derived from 5 replicates performed on 3 different days.

Highly active TNF- α and IFN- γ secretors are reproducibly detected after 3 hours of incubation in droplets

We retrieved 5 replicates of the TNF- α and IFN- γ secretion profiles from non-activated (control) and activated NK-92 cells and PMBCs. **Figure 2** shows the box plot representing the mean and standard deviation for the 5 replicates, discriminating between cells secreting TNF- α , IFN- γ , or both cytokines.

As expected, non-activated cells display a very low secretion profile, while activated cells show a high profile. The activated T cells mainly secreted TNF- α (45%) and the activated NK-92 cells IFN- γ (82%). Interestingly, almost all the T cells that secreted IFN- γ also secreted TNF- α , and almost all the NK-92 cells that secreted TNF- α also secreted IFN- γ . The standard deviation of <5% illustrates the very high reproducibility.

Multiplex and singleplex profiling results are comparable

Figure 3 shows the consistent results for the singleplex and multiplex versions of the assay.

Conclusion

The Xdrop workflow allows for cytokine secretion profiling with precise measurements at the single-cell level. The results are highly consistent and reproducible, and it is possible to investigate one or more cytokine secretion profiles simultaneously and within one working day.

Notes

Xdrop and the Xdrop DE50 Cartridge are for research use only, not for use in any diagnostic procedures.

Samplix® and Xdrop® are registered trademarks of Samplix ApS. BD Accuri™ is a trademark of Becton, Dickinson and Company. Miltenyi Biotec® is a registered trademark of Miltenyi Biotec B.V. & Co. KG. Thermo Fisher Scientific® is a registered trademark of Thermo Fisher Scientific.

Copyright © 2022 Samplix ApS

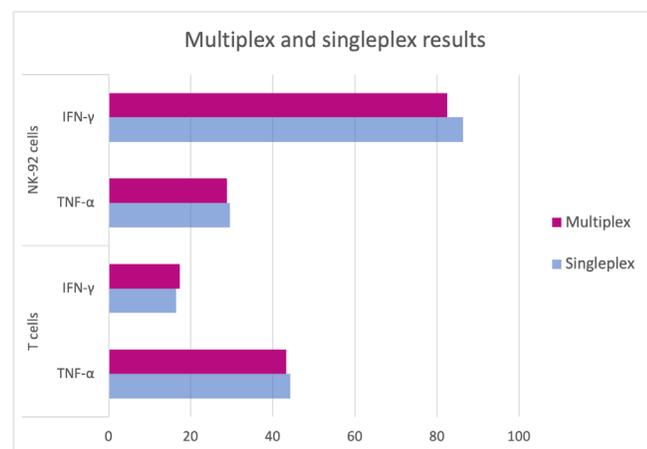


Figure 3. Clustered histogram showing the percentages of NK-92 cells and T cells secreting TNF- α and IFN- γ in multiplex and singleplex versions of the Xdrop workflow.

How Xdrop supports functional assays of mammalian cells with single-cell resolution

Using the Xdrop DE50 Cartridge, Xdrop encapsulates living mammalian cells in highly stable, ~100-picoliter, double-emulsion droplets. The picoliter reaction spaces force faster cell secretion buildup, accelerating the workflows. Xdrop processes up to 8 samples in parallel, with ~150,000 single-cell assays generated per sample in just 8 minutes. Cells within droplets can be incubated in a CO₂ incubator, analyzed on a flow cytometer, and sorted and recovered for expansion or molecular profiling.

Learn more about Xdrop at [samplix.com](https://www.samplix.com).